

# Cond 70 Vio

EC

- Resolution 0,01  $\mu$ S
- Automatic and manual temperature compensation
- Automatic calibration on 4 standard points and 1 custom
- Data logger auto and manual , up to 1000 savings in GLP
- Coloured backlit high definition LCD display
- GLP functions

## MAIN FEATURES

Professional portable meter with innovative **backlit high definition coloured LCD display**, suitable to work in any conditions, thanks also to the possibility to adjust automatically and manually contrast and brightness of the display. Possibility to select the parameters you want to show.

The **step by step guide** in setup and calibration mode allows a simple and intuitive use of the device.

**Automatic calibration of conductivity** up to 4 points and 1 point custom defined by the user.

**Adjustable reference temperature and selectable cell constant.**

**Indication with icons of the calibrated points.**

Measurement **stability indicator**.

All operations and errors are constantly monitored and shown to the operator via the **coloured LED** above the display.

**GLP function:** Date and time, recall of calibration data and function **"Calibration due"**.

**Data logger up to 1000 savings**, manual and automatic, PC connection via **micro-USB**.

**Data-Link + software** available from [www.xsinstruments.com](http://www.xsinstruments.com).

### Measuring range:

**Cond:** 0,01  $\mu$ S ...200,0 mS (0,01  $\mu$ S - automatic)


**Salinity:** 0,1 ppm ... 100 ppt

**TDS:** 0,01 mg/l ... 200,0 g/l

**Temperature:** -10...+110,0 (0,1 °C)

### Equipped with:

Conductivity standard, USB cable, power supply, user manual and carrying case.

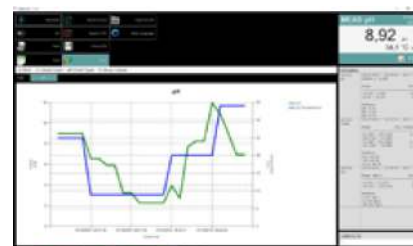

Data Link Software

Name	Description	p/n
<b>COND 70 VIO WITH CELL</b>	Complete with conductivity cell 2301T with temperature sensor integrated. Cable 1 m	50110512
<b>COND 70 VIO WITHOUT CELL</b>	With accessories, without conductivity cell	50110522